



SATURATED WIPES FOR EFFECTIVE REMOVAL OF SURFACE CONTAMINATION

Klerwipe™ 70|30 IPA Non-Sterile Pouch Wipes

Knitted, 100% polyester or low particulate, non-woven polyester/cellulose saturated with 70% v/v Isopropyl alcohol blended with Deionized (DI) Water

Klerwipe 70|30 IPA non-sterile pouch wipes are made for hard surfaces in cleanroom areas

Knitted, 100% polyester non-sterile wipes:

- ▲ Saturated with 70% v/v Isopropyl alcohol (IPA) and blended with Deionized (DI) Water
- ▲ 100% continuous filament polyester kit
- ▲ Very low particles, fibers, ions and extractables
- ▲ Basis weight: 120-140 gsm

Nonwoven, polyester-cellulose non-sterile wipes:

- ▲ Saturated with 70% v/v Isopropyl alcohol (IPA) and blended with Deionized (DI) Water
- ▲ 55% cellulose / 45% polyester hydroentangled nonwoven blend
- ▲ Low extractables and fiber and particle counts
- ▲ Basis weight: 58-78 gsm

PRODUCT QUALITY ASSURANCE

Manufactured by Berkshire in a controlled environment, Klerwipe 70|30 IPA non-sterile wipes are safe for use in core manufacturing areas.



Klerwipe™ 70|30 IPA Non-Sterile Pouch Wipes

Knitted, 100% polyester or low particulate, non-woven polyester/cellulose saturated with 70% v/v Isopropyl alcohol blended with Deionized (DI) Water

Why use Klerwipe 70|30 IPA Non-Sterile Wipes?


CONSISTENCY

Saturated Ecolab Klerwipes are absorbent which ensures a consistent application of alcohol to the entire surface area. The use of wiping can increase cleaning consistency and the quality of the final product.

OPERATOR SAFETY

Using a saturated wipe also reduces airborne chemical droplets, which reduces the risk of inhalation and improves operator safety.

Alcohol Air Concentrations After Use* Aerosol vs Trigger Spray vs Saturated Wipe

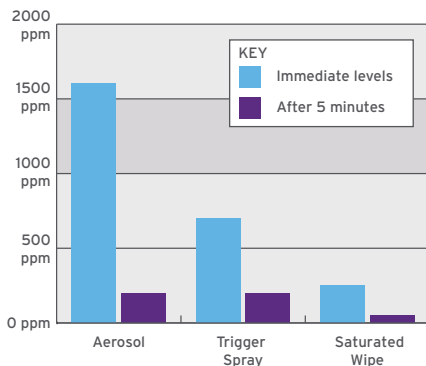

*TECHNICAL REPORT - TR9902R
Investigation into Airborne Concentrations of Alcohol when using a Trigger Spray, and Aerosol and Saturated Wipes

PRODUCT QUALITY

Klerwipes are manufactured by Berkshire in a controlled environment. These manufacturing conditions ensure Klerwipe 70|30 IPA non-sterile pouch wipes are manufactured to the highest quality and are compliant for use in your cleanroom areas.

Klerwipes are:

- ▲ Saturated with 0.2 micron filtered IPA solution
- ▲ Made with Deionized (DI) Water
- ▲ Single bagged in a resealable pouch
- ▲ Supported with Certificates of Analysis with every batch
- ▲ Supported with full documentation

PHARMACEUTICAL

Safety. Compliance. Efficiency.

APPLICATION AREA

For use on hard surfaces in cleanrooms and associated areas.

INSTRUCTIONS FOR USE

For optimum results, ensure complete wetting of surfaces. Re-seal package when not in use.

INGREDIENTS

70% v/v Isopropyl alcohol and Deionized (DI) Water.

PACKAGING INFORMATION

70|30 IPA Non-sterile pouch wipes

SUBSTRATE	UNIT SIZE	UNIT OF SALE
Knitted polyester	9 in x 9 in wipe	12 packs x 30 wipes
Woven polycellulose	9 in x 9 in wipe	12 packs x 75 wipes

FOR PRODUCT ORDER CODES AND ADDITIONAL INFORMATION, PLEASE CONTACT YOUR ACCOUNT MANAGER OR VISIT OUR WEBSITE

SHELF LIFE

Unopened, this product has a three year shelf life.

ALWAYS READ THE LABEL AND PRODUCT INFORMATION BEFORE USE.

WORLDWIDE HEADQUARTERS

1 Ecolab Place
St. Paul, MN 55102
USA

www.ecolab.com/lifesciences

